

# C.BIRD | High-throughput Screening of Recombinant Clones in Microscale

Yun Chen<sup>1</sup>, MS, Yun Chien Lin<sup>2</sup>, BS, Yu Hung<sup>2</sup>, MS, Joy Doong<sup>1</sup>, BS,  
Guan-Chiun Lee<sup>2</sup>, PhD, Charles Tsai<sup>1</sup>, PhD

<sup>1</sup>CYTENA Bioprocess Solutions, Taipei City, Taiwan

<sup>2</sup>School of Life Science, National Taiwan Normal University

## Abstract

High-throughput screening (HTS) has long been a method for selecting genetically modified clones from a pool of recombinants. To select clones with high protein and cell mass productivity for late-stage manufacturing, a need to provide comparable culture conditions between microscale and manufacture-scale has emerged. In this study, we present an innovative cell culture method for *Pichia pastoris* in 96-well plates and show a remarkable elevation in cell growth in microwells, compared to conventional static culture and aerated 96-deep-well (96-DW) shaking culture. The [C.BIRD™](#) method, which provides tunable reciprocal agitation in microwells, offers suspension culture in 96-well plates with various mixing levels and, thus, introduces higher oxygen concentration in culture medium during the cell culture process. This flexible mixing method in microscale, with no compromise on high evaporation rate, sheds light on efficient HTS methods for aerobes and overcomes current limitations on protein engineering in biopharmaceutical industries as well as academic research.

## Introduction


*Pichia pastoris*, a popular expression system for protein engineering, has found a place in commercial use in the past two decades. As a lower eukaryote, *P. pastoris* stands out because it allows for post-translational modifications as well as economical cell culture processing with high cell density outcomes. Diverse recombinant proteins are expressed in *P. pastoris*, such as antibodies, enzymes and vaccines. To attain recombinant proteins, DNA sequences of the desired proteins are integrated into the genome of *P. pastoris* for stable expression. However, this step inevitably results in clonal heterogeneity because of different numbers of expression copies integrated and off-targeted recombination. The need for an efficient method to screen for variants with high production of recombinant proteins, which are produced when *P. pastoris* reaches its stationary phase in cell growth, could not be overemphasized.

Several studies have been carried out on comparing the optimal screening methods in microscale, either in 24-well, 48-well, 96-well or 96-DW formats. But the drawbacks of each still remain a bottleneck when it comes to comparable culture conditions to those in a bioreactor, such as differences in fluidic shear rate, low oxygen transfer rate due to poor mixing level, low cell density, and high evaporation rate during the aeration process. Here, we propose the C.BIRD method, which enables tunable fluidic shear stress from choice, and, in return, promotes cell growth of aerobic organisms by increasing oxygen concentration in culture medium, with no concerns for high evaporation rates during the long-term culturing process.

## Materials and methods

The experiments were performed by Professor Guan-Chiun Lee's group at the School of Life Science of National Taiwan Normal University.

Yeast X-33 *P. pastoris* was used as the model in this study. To generate a library of stable clones, genetically modified plasmids of *P. pastoris* were integrated into the genome by electroporation. Various colonies were selected on YAG medium (20 g/L yeast extract, 5 g/L ammonium sulfate, 40 g/L glycerol) plates that contained 100 µg/mL of Zeocin, as the selection agent, and were inoculated into 1 mL of YAG medium (with 100 µg/mL of Zeocin) in 1.5 mL microcentrifuge tube. The transformants were then recovered for 3 days on a shaker at 225 rpm in a 30°C incubator.

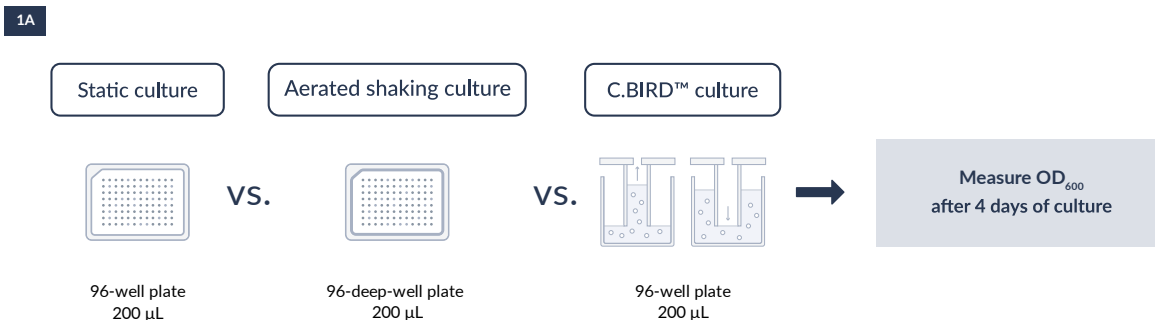


**Figure 1.** Diagram of experiment design. The measurement of OD<sub>600</sub> in each well indicated the cell growth under static, aerated shaking culture and the C.BIRD culture conditions.

For these experiments, 1 µL of each transformant was inoculated into individual wells in standard 96-well plates (or 96-DW plates), which contained 200 µL of YAG medium (with 100 µg/mL of Zeocin) and were incubated at 30°C under either static culture, aerated 96-DW shaking culture at 225 rpm, or C.BIRD culture with 13 seconds of mixing period. After 4 days of culture, absorbance at wavelength of 600 nm was measured in 96-well plates by Thermo Scientific's Multiskan GO Microplate Reader.

## Results and discussion

The experiments were designed as shown in **Figure 1**, and the absorbance in each well of static culture, aerated 96-DW shaking culture, and C.BIRD culture was measured by microplate reader on day 4. The results showed that the OD<sub>600</sub> in the static group remained low at an average of 0.14, while that of the C.BIRD and aerated 96-DW shaking groups were at 0.43 and 0.29, respectively (**Figure 2A-C**). In the same clones, the C.BIRD group got 3.10-fold and 1.55-fold augmentation in OD<sub>600</sub> value on average, compared to static and aerated 96-DW shaking groups. This indicated that the C.BIRD cell culture method could efficiently promote the cell growth of *P. pastoris*, compared to other conventional culture methods. Furthermore, the C.BIRD culture method provided a way to select the potential high-growth rate clones in microscale, such as clones A7, A8 and B3 in **Figure 2C**, and enabled an efficient and low-labor HTS method with more comprehensive outcomes. Additionally, with the overall elevation of OD<sub>600</sub> value in the C.BIRD group, more clones that


2

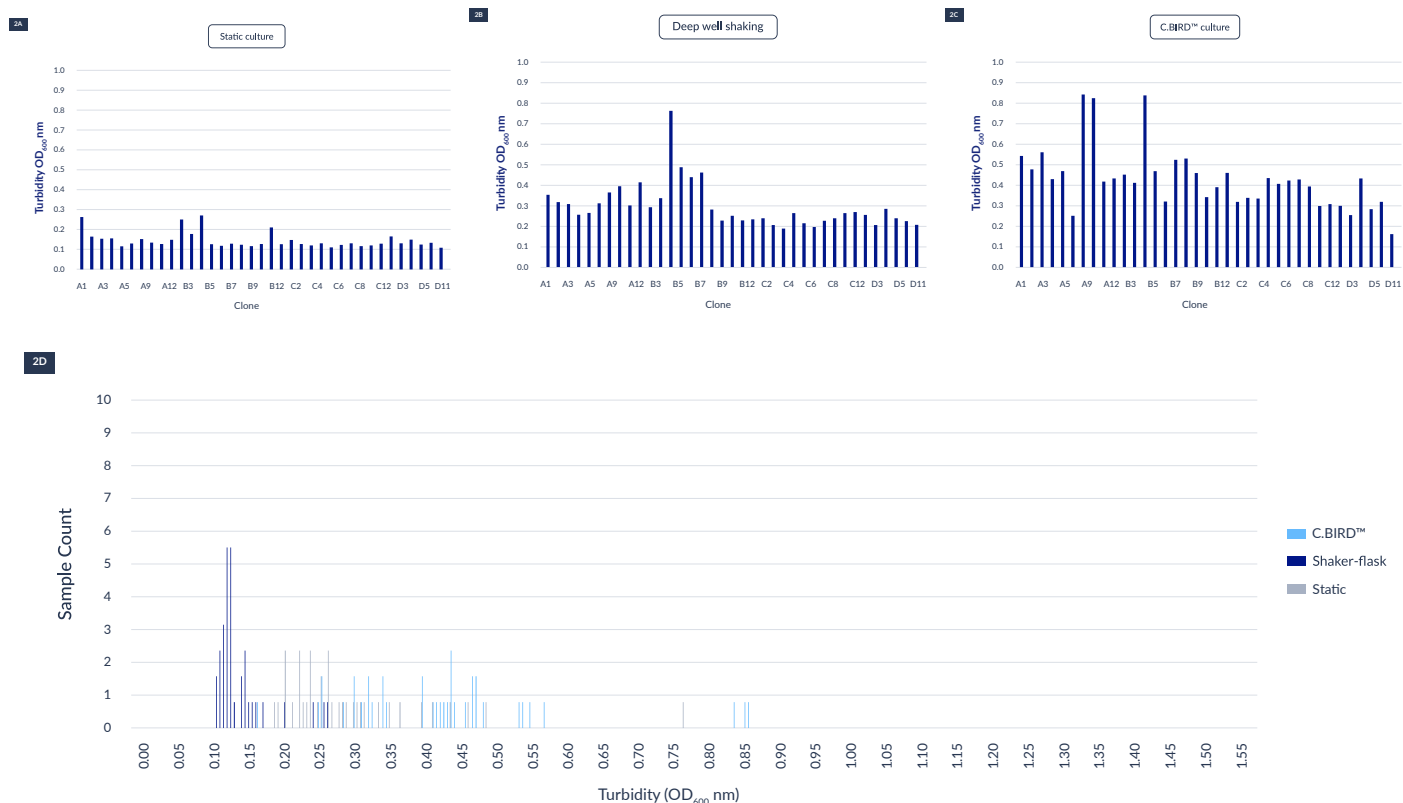

**Figure 2.** The cell densities under three culture methods after 4 days of culture. The x-axis represented the name of the clones. **A)** static culture, **B)** aerated 96-DW shaking culture, and **C)** C.BIRD culture were shown. **D)** Integrated data showed the distribution of turbidity and its sample counts of the three culture methods.

harbored high cell mass productivity were observed, compared to other culture method groups. This showed that the C.BIRD culture method could not only significantly increase cell growth in *P. pastoris* but also help selecting potential highly productive clones from the pool.

In parallel, a cross-contamination test was performed to secure the independence of each clone as well as the following measurements. Different clones were inoculated into 200  $\mu$ L of YAG culture medium in the wells of a 96-well plate, while 200  $\mu$ L of sterile culture medium was placed in the adjacent wells without inoculum, as shown in **Figure 3A**. After 4 days of continuous culture with the C.BIRD device, it showed that the OD<sub>600</sub> value of the sterile wells remained the same as on day 1, which indicated that there was no well-well cross-contamination (**Figure 3B**). On the other hand, the wells, inoculated with transformants, showed increasing OD<sub>600</sub> values during 4 days of culture, therefore demonstrating the capability for cell division of the clones (**Figure 3C**).

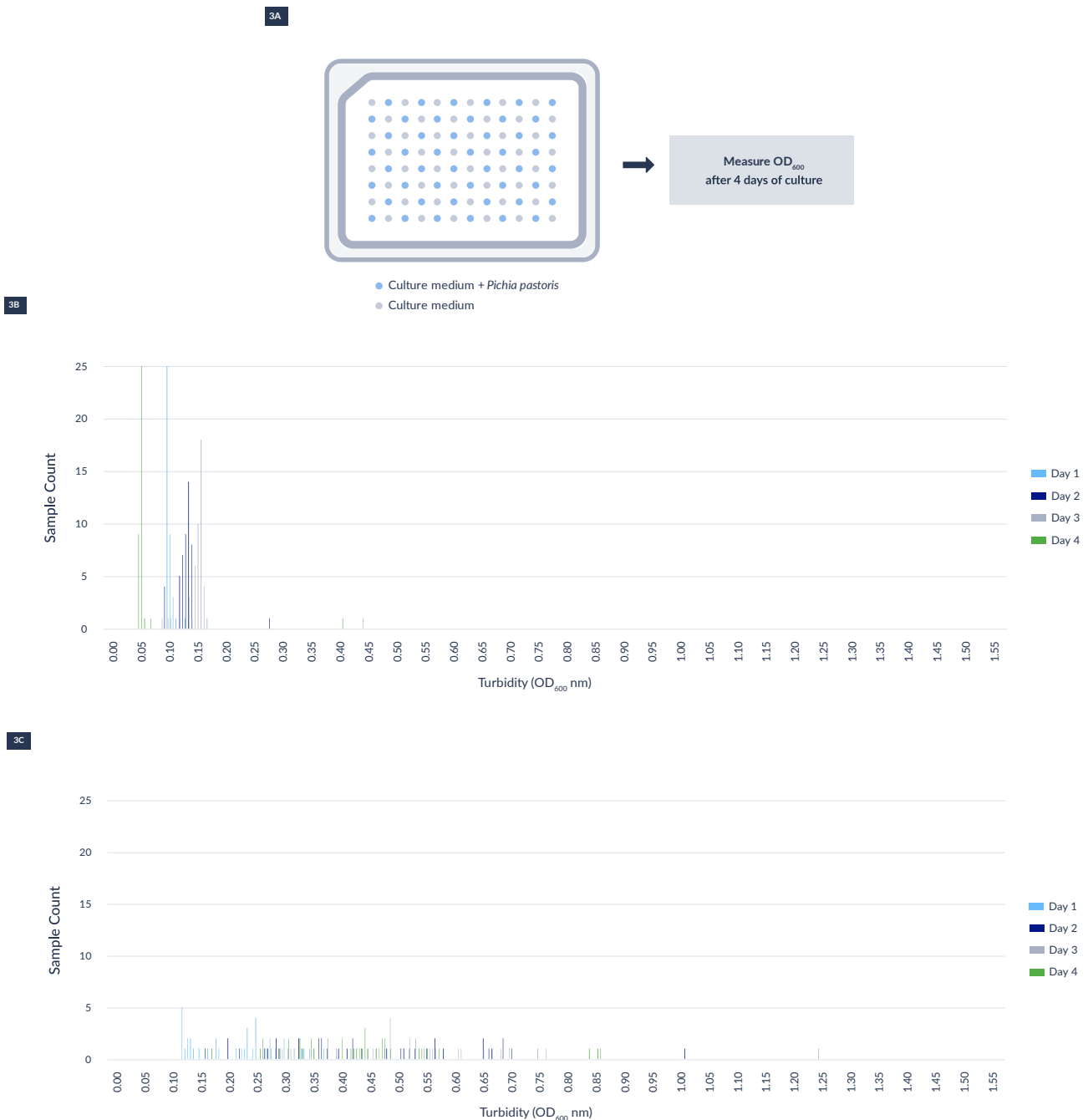

**Figure 3.** Cross-contamination test results. **A)** The experiment design of cross-contamination test. **B)** In culture medium group, it showed that there was no cross-contamination between adjacent wells, while **C)** the cell division of *P. pastoris* remained active.

## Conclusion

Introducing the C.BIRD method to standard 96-well plates provided higher oxygen transfer rates, which furnished the aerobic cell culture with better culture conditions. Our results show that efficient mixing during the *P. pastoris* cell culture process dramatically and evenly elevated cell density, allowing researchers to select highly productive clones in a high-throughput fashion, with no need for laborious process. Moreover, the results also ensured the independence of each well during long-term culture. This undoubtedly fortified the strength and reliability of the C.BIRD method on selecting clones from genetically modified transformants.

The C.BIRD method refined the conventional HTS methods in *P. pastoris* by enabling a well-mixed and aerated culture condition in 96-well plates and increased the average cell density of each clone in comparison with other methods. The C.BIRD method provides an ideal cell culture condition for aerobic organisms and offered a more efficient way for HTS in *P. pastoris*.

## Reference

Kaushik N, Lamminmäki U, Khanna N, Batra G. Enhanced cell density cultivation and rapid expression-screening of recombinant *Pichia pastoris* clones in microscale. *Scientific Reports*. 2020; 10(1): 7458. [DOI:10.1038/s41598-020-63995-5](https://doi.org/10.1038/s41598-020-63995-5).

## Contact us

Phone: +886-2 2720 6135

Email: [info@cytena-bps.com](mailto:info@cytena-bps.com)

Website: [www.cytenua-bps.com](http://www.cytenua-bps.com)



© 2021 CELLINK AB. All rights reserved. Duplication and/or reproduction of all or any portion of this document without the express written consent of CELLINK is strictly forbidden. Nothing contained herein shall constitute any warranty, express or implied, as to the performance of any products described herein. Any and all warranties applicable to any products are set forth in the applicable terms and conditions of sale accompanying the purchase of such product. CELLINK provides no warranty and hereby disclaims any and all warranties as to the use of any third-party products or protocols described herein. The use of products described herein is subject to certain restrictions as set forth in the applicable terms and conditions of sale accompanying the purchase of such product. CELLINK may refer to the products or services offered by other companies by their brand name or company name solely for clarity and does not claim any rights to those third-party marks or names. CELLINK products may be covered by one or more patents. The use of products described herein is subject to CELLINK's terms and conditions of sale and such other terms that have been agreed to in writing between CELLINK and user. All products and services described herein are intended FOR RESEARCH USE ONLY and NOT FOR USE IN DIAGNOSTIC PROCEDURES.

The use of CELLINK products in practicing the methods set forth herein has not been validated by CELLINK, and such nonvalidated use is NOT COVERED BY CELLINK'S STANDARD WARRANTY, AND CELLINK HEREBY DISCLAIMS ANY AND ALL WARRANTIES FOR SUCH USE. Nothing in this document should be construed as altering, waiving or amending in any manner CELLINK's terms and conditions of sale for the instruments, consumables or software mentioned, including without limitation such terms and conditions relating to certain use restrictions, limited license, warranty and limitation of liability, and nothing in this document shall be deemed to be Documentation, as that term is set forth in such terms and conditions of sale. Nothing in this document shall be construed as any representation by CELLINK that it currently or will at any time in the future offer or in any way support any application set forth herein.